
# T1-173 & T1-174 Transmitter (AA)

## Instructions For Use & Care

English  
Inglés

### OPERATING INSTRUCTIONS

1. Install two (2) AA batteries. If you're using rechargeable batteries, they must be charged before using.
2. Plug the microphone cord into the "Mic" jack on top of the Transmitter.
3. Place the Transmitter in the belt clip case provided.

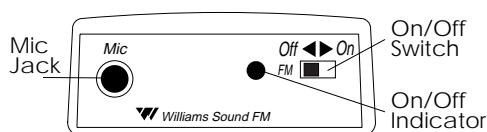4. Turn the power switch on top of the Transmitter to "On."
5. The microphone should be placed as close to the speaker's mouth as is practical. For lapel mics, attach the microphone to a collar, lapel, or tie.

The transmitter can be placed in a pants pocket, or clipped onto a belt or waistband. Make sure the Transmitter is turned OFF when not in use.

#### *Note On The Transmitter Antenna:*

The microphone cord is the transmitting antenna. Do not bunch up the cord or wrap it around the transmitter. For maximum range, the cord should hang as straight as possible.

### BATTERY INSTRUCTIONS


This transmitter can use disposable AA batteries (alkaline or lithium) or a Williams Sound rechargeable battery (BAT 026). We recommend using BAT 026 AA Ni-Cad or BAT 001 AA Alkaline batteries. BAT 026 batteries will last about 8 hours per charge when used with Williams Sound Transmitters. BAT 001 batteries will last approx. 12 hours when used with Williams Sound Transmitters. Batteries from other suppliers may provide different operating life.

### EXPECTED BATTERY LIFE

Alkaline batteries (BAT 001):	12 hrs
Rechargeable Ni-Cad Batteries (BAT 026):	8 hrs per charge

### INSTALLING AND REPLACING BATTERIES

1. Open the battery compartment by lifting the tab on the bottom of the transmitter with a coin or finger.
2. To remove depleted batteries, pull up on the fabric strip.

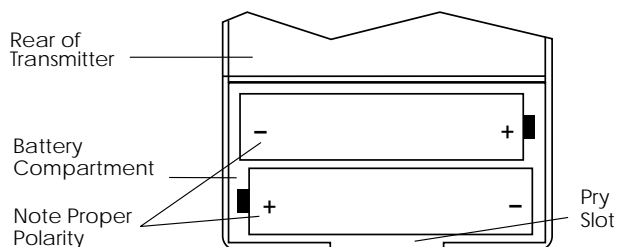3. To install new batteries, press them into place over the fabric strip. Be sure to observe proper polarity (+/-).
4. Close the lid.
5. When the sound becomes weak or distorted, replace or recharge batteries. If recharging, use the Williams Sound CHG 200 or CHG 1600 Battery Charger.

### WARRANTY NOTICE

Your Williams Sound Transmitter is engineered and designed to provide you with many years of reliable service. Williams Sound warrants it against defects in materials and workmanship for FIVE (5) years EXCEPT FOR microphones, batteries, chargers, cables, carry cases, and all other accessory products. Accessory products carry a 90 day warranty.

If the transmitter fails within the warranty period, Williams Sound will determine whether to repair or replace it. This warranty does not apply to physical damage, transmitters that have been abused or mis-used, or modified.

Williams Sound Products are NOT designed for use in extreme temperature, humidity, or chemical environments. The introduction of chemicals such as chlorine, salt water, or human sweat into the product will cause damage to the circuitry. Damage due to these causes IS NOT covered under the product warranty.

 **Williams Sound**<sup>®</sup>  
Helping People Hear  
[www.tour-guide-systems.com](http://www.tour-guide-systems.com)

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